




Operators manual

NT-851ABR – Adjustable backrest



NSN: 6530-25-160-2915

Before use on the NODIN Medevac Stretcher (NT-820MVU)

Remove the shoulder straps from the two webbing-loops between the topmost hinges on the stretcher


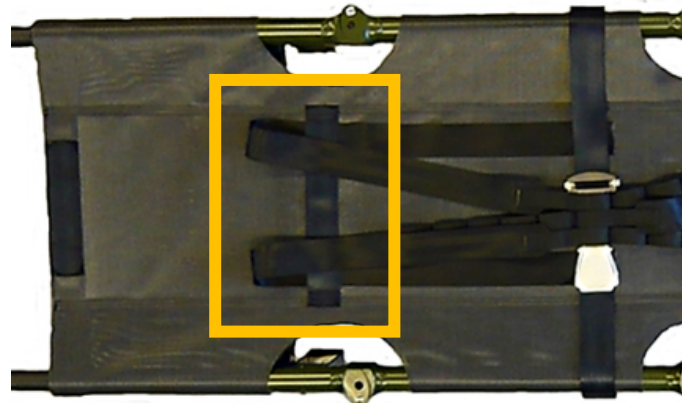




Attaching the backrest to the NODIN Medevac Stretcher

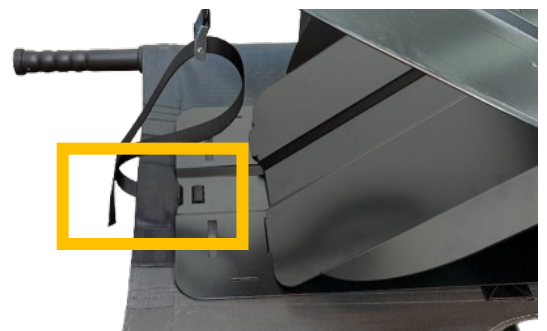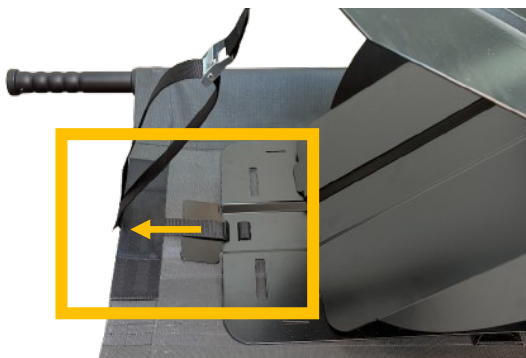

Place the Backrest on the stretcher and push the web strap **through** the loop in the stretcher head end and pull the flap on the backrest bottom plate into the loop for secure fastening.

Insert the shoulder straps into the two slots at the lower end of the backrest, closest to the shoulder strap anchoring point.

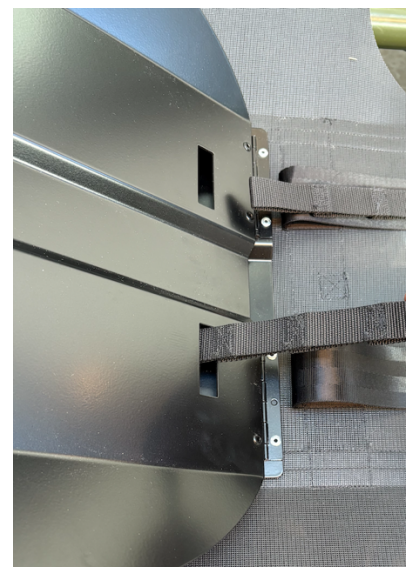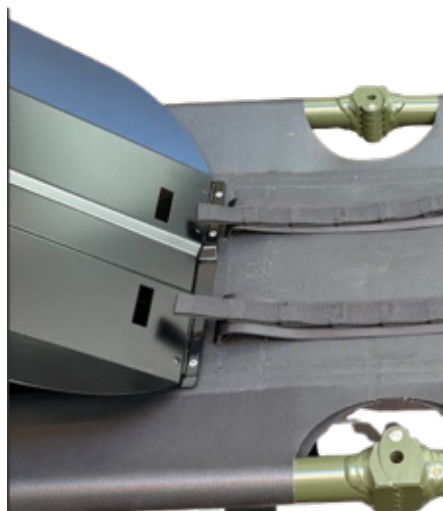

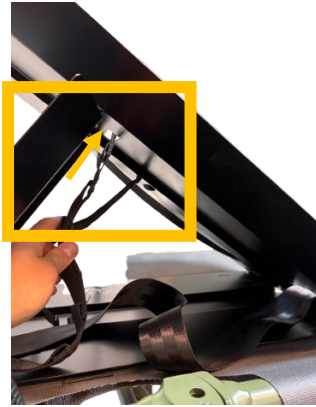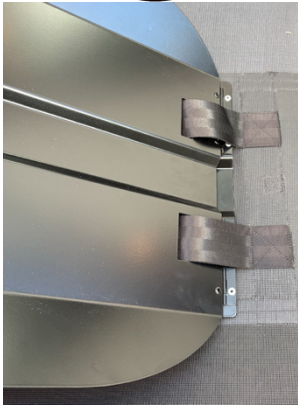

Feed the straps through and along the rear of the backrest.

Insert the shoulder straps through the slots at shoulder level on the backrest and feed them forwards towards the upper body buckle

Make sure to adjust for any slack when adjusting the angle of the backrest.

Adjusting the angles



Adjust the back-plate support to the desired angle and tighten the web strap slightly to fix the position/angle.

Loosen the web strap to adjust the angle.



NOTE:

Be aware that if the patient is securely strapped into the harness on the NODIN Stretcher, the straps are fastened to the backrest as described in this manual AND the webstrap at the top end is tightened, the patient will NOT be able to sit up without help.