


# Springtime™

## Self-ligating bracket

- Interactive Self-ligating bracket with a clip
- Made of a hard 17-4 stainless steel material.
- Metal injection molding (MIM) advanced technology
- Contoured appearance for maximum comfort
- Mesh base with laser mark
- Available in Roth .022.018 and MBT will be available next year

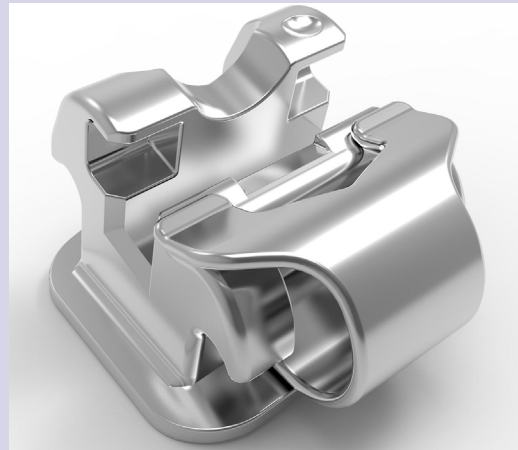

### Passive Phase

Small round wires slide freely, initiating the tooth movement process as the archwire gently levels the teeth and coaxes them into alignment.



### Expressive Phase

Square or rectangular wires are seated into the base of the slot without contacting the clip. Rotations are corrected and space closures are completed.



### Active Phase

Rectangular archwires extend beyond the slot to fully engage the clip, providing the active control necessary for finishing uprighting of the roots and adjusting the torque

**GESTENCO**  
*Smiles forever*

Gestenco International AB, P.O. Box 240 67, SE-400 22 Gothenburg, Sweden  
Phone +46 (0)31 81 00 35. Fax +46 (0)31 81 46 55, [www.gestenco.com](http://www.gestenco.com)

## Springtime™ Formula Roth\*

Type	Description	Torque	Ang.	Roth .022
Centrals	Upper right	+12	5	SP-11-55
	left			SP-21-55
Laterals	Upper right	+8	9	SP-12-55
	left			SP-22-55
Cuspids w/hook	Upper right	-2	10	SP-13-55HK
	left			SP-23-55HK
1st Bicuspid	Upper right	-7	0	SP-14-55
	left			SP-24-55
1st Bicuspid	Upper right	-7	0	SP-14-55HK
	left			SP-24-55HK
2nd Bicuspid	Upper right	-7	0	SP-15-55
	left			SP-25-55
2nd Bicuspid	Upper right	-7	0	SP-15-55HK
	left			SP-25-55HK
Anteriors	Lower right	-1	0	SP-40-55
	left			SP-30-55
Cuspids w/hook	Lower right	-11	7	SP-43-55HK
	left			SP-33-55HK
1st Bicuspid	Lower right	-17	0	SP-44-55
	left			SP-34-55
1st Bicuspid	Lower right	-17	0	SP-44-55HK
	left			SP-34-55HK
2nd Bicuspid	Lower right	-22	0	SP-45-55
	left			SP-35-55
2nd Bicuspid	Lower right	-22	0	SP-45-55HK
	left			SP-35-55HK

**GESTENCO**  
*Smiles forever*

Gestenco International AB, P.O. Box 240 67, SE-400 22 Gothenburg, Sweden  
Phone +46 (0)31 81 00 35. Fax +46 (0)31 81 46 55, [www.gestenco.com](http://www.gestenco.com)