

iFlash Test Panels

Main Stream Hormonology	Reproductive Health	Autoimmunity	Infection
Thyroid TSH Anti-TPO Anti-TSHR Tg Anti-Tg FT3 FT4 T3 T4 Tumor Markers PG I PG II AFP CEA CA 125 CA 15-3 CA 19-9 Free PSA Total PSA CYFRA 21-1 NSE Anemia Ferritin Folate Vitamin B12 Cardiac Markers CK-MB Myoglobin Troponin-I BNP D-Dimer ** Bone Metabolism 25-OH Vitamin D PTH Metabolism Insulin C-Peptide Cortisol ACTH Adiponectin Liver Fibrosis HA PIIIPN-P CoI IV Laminin	Fertility AMH Inhibin B Inhibin A LH FSH HCG Prolactin DHEA-S Unconjugated Estriol E2 PAPP-A free βHCG Progesterone Testosterone ToRCH Toxo IgG Toxo IgM CMV IgG CMV IgM HSV-1 IgG HSV-2 IgG HSV-1 IgM HSV-2 IgM Rubella IgG Rubella IgM Anti-Phospholipid Syndrome Anti-Cardiolipin Cardiolipin IgM Cardiolipin IgG Cardiolipin IgA Anti-β2-Glycoprotein I β 2-Glycoprotein I IgM β 2-Glycoprotein I IgG β 2-Glycoprotein I IgA * Without CE ** Coming Soon	Rheumatoid Arthritis Anti-CCP RF IgG RF IgM RF RA33 IgG Connective Tissue Disease ANA dsDNA IgG Sm IgG SS-A IgG SS-B IgG Scl-70 IgG Jo-1 IgG RNP70 IgG Autoimmune Liver SMA IgG AMA-M2 Vasculitis Diagnostics PR3 IgG GBM IgG MPO IgG Autoimmune Diabetes IA-2A IAA GADA ICA ZnT8A Celiac Disease tTg IgA ** β 2-Glycoprotein I IgM tTg IgG ** DGP IgA ** DGP IgG **	EBV EB VCA IgG EB VCA IgM EB VCA IgA RF EB EA IgM EB NA IgG EB NA IgA Respiratory Disease Mycoplasma pn. IgG Mycoplasma pn. IgM Chlamydia pn. IgG Chlamydia pn. IgM Blood Bank Anti-TP (Syphilis) Anti-HCV * HBsAg * Anti-HBs * HBeAg * Anti-HBe * Anti-HBc * HIV Combo * Inflammation PCT IL-6 ** SARS-CoV-2 2019-nCoV IgM 2019-nCoV IgG 2019-nCoV IgA * 2019-nCoV Ag 2019-nCoV NAb SARS CoV-2 IgG-S



Specifications


General Information	
Throughput	Up to 180 tests/hour
Fastest time to first result	12 minutes
Detection principle	Magnetic particle-based direct chemiluminescence
Chemiluminescence label	Acridinium Ester
Sample Supply System	
Sample position	5 samples per rack, 50 sample positions with STAT function
Sample volume	5~150 μL
Reaction vessels loader capacity	1000 reaction vessels reserve for continuous loading
Reagent System	
Reagent position	20 positions with onboard refrigeration
Reagent carousel temperature	2-8°C
Reagent detection	Reagent level tracking
Reagent pack	2x50 tests per kit
Auto dilution	Maximum of 1:400
Reaction Carousel	
Capacity	75 reaction vessels
Incubation temperature	37.0°C±0.3°C, with a fluctuation of 0.1°C
Mixing system	Non-contact eccentric vortex mixer
Magnetic Separation System	
Washing method	4-step magnetic separation
Temperature	37.0°C±0.3°C, with a fluctuation of 0.1°C
Operating Conditions	
Power supply	100~240 V, 50/60 Hz
Power consumption	1000 VA
Temperature	10~30°C
Ambient humidity	20%-70% RH
Physical dimensions	1030 mm (L) × 750 mm (W) × 710 mm (H)
Weight	150 kg
Consumables	
Product No.	C89997/C80001/C89949/C89948/C6100
Description of Goods	Wash Buffer (1x), Wash Buffer (10x), Pre-Trigger Solution (1x) Trigger Solution (1x), Reaction Vessels
Spec.	10L, 1Lx4 bottles, 220mLx4bottles, 220mLx4bottles 2000 cuvettes / box



SHENZHEN YHLO BIOTECH CO., LTD.

Headquarter: Building 1, YHLO Biopark, Baolong 2nd Road, Baolong Subdistrict, Longgang District, 518116 Shenzhen, P.R.China

Tel: +86 755 26601910
 Email: marketing@szyhlo.com
 Website: www.szyhlo.com

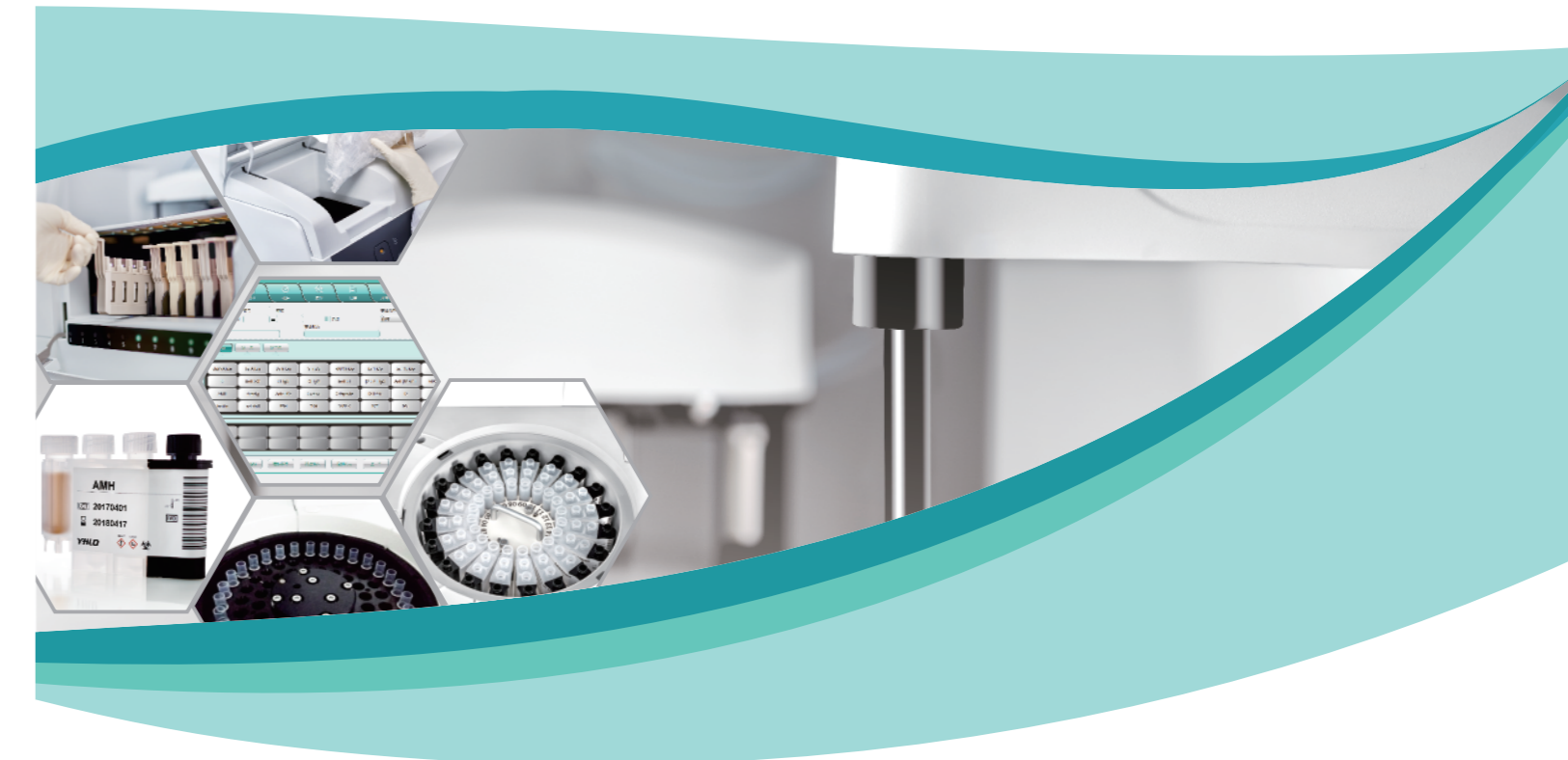

iFlash 1800
 Chemiluminescence
 Immunoassay Analyzer